



Medit Intraoral Scanner Lineup

MEDIT



Digital Dentistry, *Your Way*

Maximize efficiency and precision with Medit

Infinite Possibility

At Medit, we are not just about creating products; we are re-imagining possibilities.

Our commitment to innovation and customer care propels the digital dentistry revolution, offering 3D digital solutions that transform dental clinics and labs.

With our award-winning intraoral scanners, based on our own patented state-of-the-art technology and collaborative software, we are setting new standards in dental technology.

Our Medit Community

With continuous feedback, active discussion, and helpful suggestions, the Medit community is always connected as customers and colleagues actively support and help one another. There's a reason why over 100,000 dental professionals around the globe trust us.

Scan the QR code to join the Medit revolution worldwide and experience it yourself!



User Group



Medit Academy



Help Center





MEDIT *i900* classic
MEDIT *i900*

Effortless Precision,
Your Way



Effortless Precision, *Your Way*

Advanced scanning technology available in both
Touch and Button interfaces.



Fast, Easy Scanning for an Effortless Workflow

- Capture a larger field of view in one go, enabling quicker scans.
- Easily scan edentulous patients for efficient treatment.



Compact, Lightweight and Comfortable to Handle

- Weighing just 165g, lighter than a smartphone to reduce wrist strain during extended use.
- Slim, ergonomic design for a comfortable and secure grip.



Seamless Hygienic Design That's Easy to Clean

- Seamless, unibody design for easy cleaning.
- Built-in UV sterilization keeps the device clean from the inside out.



Tactile Button

Click and scan with
confident precision

Wand Tip

Multiple tip sizes to
customize patient comfort

Intelligent Sensor

LED and haptic
feedback

Light & Compact

165g wand weight
223mm compact size

3rd-Generation Optical Engine

Advanced 10-bit
camera technology

LED Indicator

Stay Informed from
Connection to Scanning

Fast Preheating

Ready to scan
instantly

Touch Band

Real-time, light and
haptic feedback

Touch Pad

Customizable command actions
Mouse-free with smart interaction

Unibody design

Jointless design with
perfect balance
No biohazard concerns
Reinforced, hassle-free
cable

Medit i700: The Lightweight Powerhouse

Fast, precise scanning - making visits quicker and more comfortable without sacrificing accuracy. With both powerful hardware and intelligent software, the Medit i700 is the key to unleashing your clinic's full potential.



180° Reversible Tips

Remote Control Mode

Detachable Cable

245g Featherweight

UV-C LED Disinfection

2X Super Fast





No More Wires

Free from wires, scanning is comfortable at any angle. Scan proximal areas effortlessly without constraints in movement.



Simply Switch & Scan

Scan effortlessly across multiple chairs. With additional hub setup, easily scan across multiple chairs with a single click.



Supercharged Yet Light

Outstanding weight balance evenw tih a battery attached for added comfort and uninterrupted operation with long battery life (three batteries included).




Powerful Hardware Meets Versatile Software

Optimized for convenience. Not only from hardware aspects, but also software. *i700* wireless, packed with the latest of MEDIT technology.

Beyond *i700*'s Proven Tech

Still the same **super fast**, **light**, and **accurate *i700***. But now **wireless**.

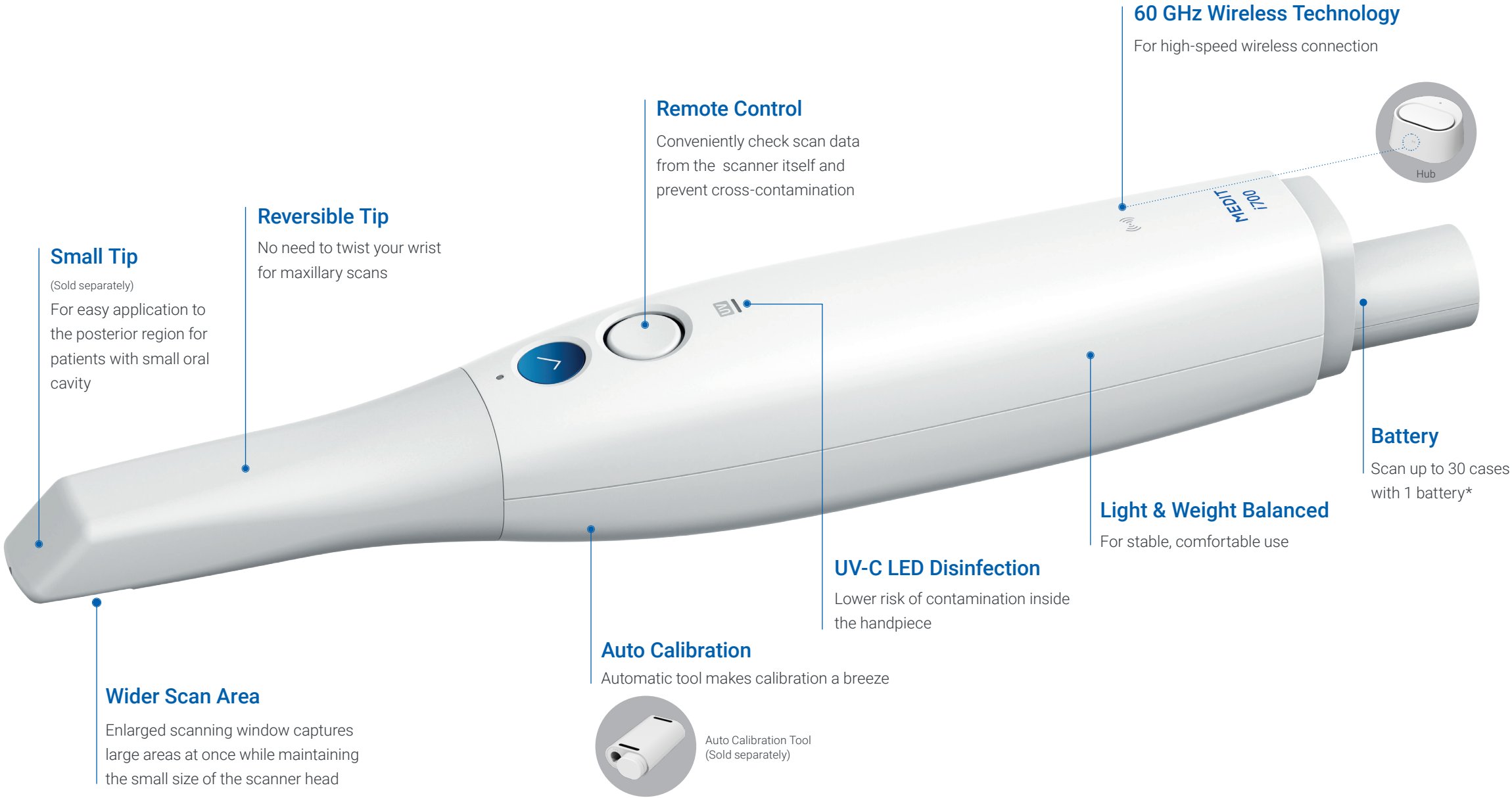

60 GHz Wireless Technology

For high-speed wireless connection

Remote Control

Conveniently check scan data from the scanner itself and prevent cross-contamination

Reversible Tip

No need to twist your wrist for maxillary scans

Small Tip

(Sold separately)

For easy application to the posterior region for patients with small oral cavity

Wider Scan Area

Enlarged scanning window captures large areas at once while maintaining the small size of the scanner head

Auto Calibration

Automatic tool makes calibration a breeze



Auto Calibration Tool (Sold separately)

UV-C LED Disinfection

Lower risk of contamination inside the handpiece

Light & Weight Balanced

For stable, comfortable use

Battery

Scan up to 30 cases with 1 battery*

When the scanner's cable disappears, so do restrictions to scan movements. Strengthened impact resistance minimizes damage from accidental drops. Weight balanced design of the *i700* wireless provides a stable, comfortable scanning experience.

*Performance may vary according to usage conditions.

The Brain of the *i700* with a Splash of Color

**Super Accurate,
Super Light,
Super Smooth.**

The digital workflow begins with accurate data.
Medit *i600* precisely acquires data for effortless digital dentistry.

***High accuracy**

Full-arch accuracy: $10.9\ \mu\text{m} \pm 0.98$

*Acquires accurate scan data from anterior segment to the occlusal surface

Astonishing Scanning Experience

Light,

245g featherweight with improved grip for comfort.

Smooth,

Improved speed and optimal scanned images deliver a smoother scanning experience.

Vivid,

Transmits vivid scan data via the new optical system used in the Medit *i700*.

Scan Swiftly with Confidence

Take full control of your scanning with clicks of one button.
No need for numerous cables and hubs!
Simply Plug & Scan using a single cable for direct connection with a PC.



Medit Link

Effortless clinic data management with exceptional results.



Create and manage scan data matched with patient information in the electronic chart program and link them with treatment solutions of various workflows. The various apps within Medit Link help provide a better consultative experience that is easier for patients to understand.


- Encrypted cloud storage for reliable protection against cyber security attacks.
- Enabled data sharing and real-time communication between dentists, patients, and lab technicians.
- Accessible anytime, anywhere on a PC, tablet, or mobile device with internet connection.

Over the years, Medit has tirelessly refined its software, culminating in a comprehensive suite of industry-leading digital dentistry applications. Offering unparalleled value to our users, Medit has set the standard for innovation and functionality, earning recognition for the best IOS software applications in the field.

Download **Medit Link**



Digital Dentistry with **Medit Link** & **Medit Apps**

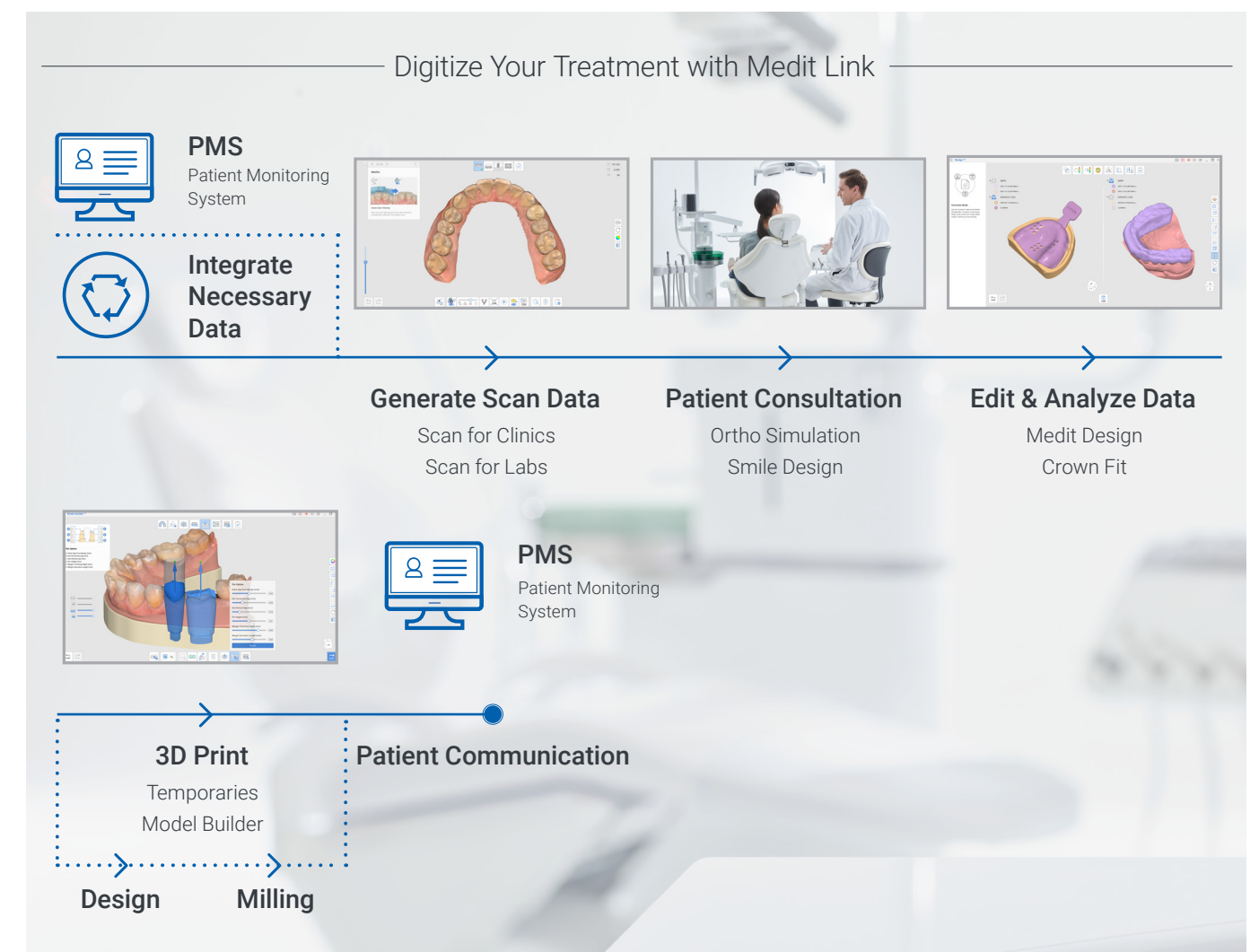

Use Your Scan Data for Much More

Scan data created with the Medit scanners can be linked to the workflow of various solution providers via Medit Link. The data can also be used with any software that can open standard 3D formats such as STL, PLY, and OBJ as we continue to expand the utilization of our open system.

What you can do in Medit Link :

Clear aligners | Crown and bridges | Implant-supported restorations | Full and partial dentures | Orthodontics | Implant surgical guides | Sleep apnea devices | Dental models | 3D printing and more.








Medit Apps

Medit Apps quickly introduce digital workflows using scan data and help you utilize them in a diverse range of cases. You can analyze data, create restorations and models, and communicate with your patients in the simplest way. With Medit Apps, a lot more can be done on the chairside like designing splints, temporaries, bridges, and so on! Additionally, there is no subscription cost or annual fee to be a part of the Medit ecosystem.

Medit Intraoral Scanner Lineup Comparison

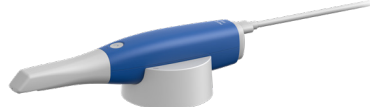

Category		i900 classic		i900	i700 wireless	i700	i600
Model							
Scanning Technology	Imaging Technology	3D-in-motion video technology / 3D full color streaming capture					
	Anti Fogging Technology	Adaptive anti-fogging					
Accuracy	Full Arch	10.9 μm ± 0.98*					
Handpiece	Interface Type	Button		Touch		Button	
	Dimensions	223 x 37 x 35 mm			313 x 44 x 47 mm	248 x 44 x 47 mm	
	Total Weight	165g			328g	245g	
Tip	Scan Area	18 x 15 mm			15 x 13 mm		
	Autoclavable	Sterilize for 10 minutes at 135°C (275°F) at gravity type and dry for 30 minutes. Sterilize for 4 minutes at 134°C (273°F) at pre-vacuum type and dry for 20 minutes.					
	Reversable Tip	No			Yes		Yes
Special Mode	Remote Control Mode	No		Yes			No
	UV-LED Disinfection	Yes					No
	Smart Shade Guide	Yes					
	HD Mode	Yes					No
MAC OS Compatibility		Yes					
Connectivity	Connection Type	Wired USB 3.0					
	Plug and Scan	Yes					
	Detachable Cable	Yes					

* May vary depending on scan object

